



GAIT BELTS

	SINGLE PATIENT USE
	NON-STERILE
	NOT MADE WITH NATURAL RUBBER LATEX
RX ONLY	FEDERAL U.S.A. LAW RESTRICTS THIS DEVICE TO SALE OR USE BY OR ON THE ORDER OF A PHYSICIAN OR PROPERLY LICENSED PRACTITIONER.

IMPORTANT INFORMATION

Please read all instructions, warnings, and cautions before use. Correct application is essential for proper functioning of the product. Use only on the person it was provided to by a healthcare professional and only for the use it was intended.

INTENDED USE

The DeRoyal® Gait Belts are intended to provide caregivers a secure point to hold while assisting patients with walking and/or transfer.

CONTRAINDICATIONS


Never use on a patient if an I.V. or wound site could be compromised by the device.

! WARNINGS

- Discontinue use and consult your prescriber if the product or its components break or become damaged, or cleaning is no longer sufficient to sanitize product.
- Inspect the fit of the product on the patient to ensure fit is proper. **DO NOT OVERTIGHTEN.** Check regularly to ensure breathing and circulation are not compromised.

CAUTIONS

- A physician or properly licensed practitioner (a "prescriber") who is familiar with the use and purpose of this product must fit it to the user. The prescriber has a duty to provide wearing instructions and risks related to the use of this product to other healthcare practitioners treating the users and the users themselves, including duration of use. The instructions provided in this sheet do not supersede hospital protocol or direct orders of the prescriber.
- Belt should be worn over patient's clothing. Consult prescriber immediately if patient experiences sensation changes, unusual reactions, swelling or increased pain while using this product. Discontinue use if


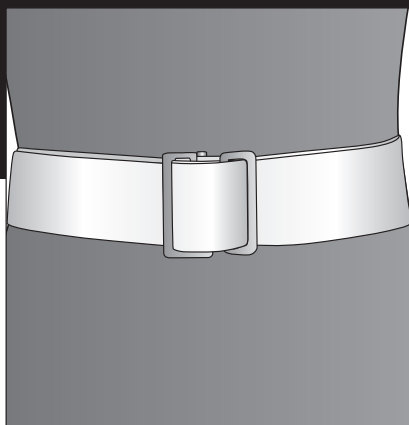


pressure injuries develop.

- Take special care if the user is diabetic or has poor circulation. These users may have decreased skin sensitivity and are at greater risk for poor peripheral circulation and pressure injuries.



INSTRUCTIONS FOR USE

1. Measure the patient's waist circumference and select appropriate size.
2. Wrap the belt around the patient's waist with the buckle closure in front. **NOTE:** Ensure the label is on the outside of the device.
3. Secure the buckle closure. See below (3A or 3B) for buckle application instructions. **NOTE:** Do not over tighten. Two fingers should fit between the device and the patient.
4. Once the closure is secure, buckle should be positioned to the back of the patient so that the patient cannot release during transfer. Tuck excess belt out of patient's reach.

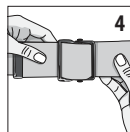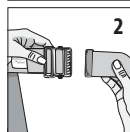

TIGER-TOOTH BUCKLE CLOSURE

3A. Feed the metal tip of the belt through the portion of the buckle with the tiger-teeth. Slide belt over the tiger-teeth, and feed the end of the belt through the empty slot on the opposite side of the buckle. Pull to tighten. Make sure buckle and belt lay flat against the patient.

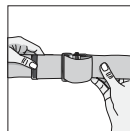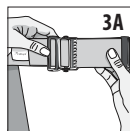

QUICK-RELEASE BUCKLE CLOSURE

3B. Insert the male end into the female end of the buckle. Pull to tighten.



SIZING

DIMENSIONS

54"L x 2"W (137cm x 5cm)

56"L x 2"W (142cm x 5cm)

60"L x 2"W (152cm x 5cm)

70"L x 2"W (178cm x 5cm)

84"L x 2"W (213cm x 5cm)

CLEANING AND/OR MAINTENANCE

COTTON GAIT BELTS: Machine wash with mild fabric detergent and warm water. Tumble dry on low. Do not bleach. **NOTE:** Fasten buckles before washing. Ensure buckles are working before re-using the product. Product may experience some shrinkage from multiple washes. If the product is damaged or cleaning is insufficient, do not use it and replace it with a new one.

STORAGE AND TRANSPORT CONDITIONS

	KEEP DRY
	KEEP AWAY FROM SUNLIGHT

In addition to the competent authority in the country where the patient resides, serious incidents must be reported to DeRoyal Industries, Inc.

WARRANTY

DEROYAL® PRODUCTS ARE WARRANTED FOR NINETY (90) DAYS FROM THE DATE OF SHIPMENT FROM DEROYAL AS TO PRODUCT QUALITY AND WORKMANSHIP. **DEROYAL'S WRITTEN WARRANTIES ARE GIVEN IN LIEU OF ANY IMPLIED WARRANTIES, INCLUDING WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE.**



Manufactured by:
DeRoyal Industries, Inc.
200 DeBusk Lane
Powell, TN 37849 USA
888.938.7828 or (001) 865.938.7828
www.deroyal.com

PART#0-700 | REVISED 5/2023

©2023 DEROYAL INDUSTRIES, INC.

All Rights Reserved. DeRoyal and the DeRoyal logo, are trademarks or registered trademarks of DeRoyal Industries, Inc.